



Your Trio **Siltac**[®] and Trio **Silvex**[®] sample guide

—
Feel the difference, experience the benefit
of our Responsive Silicone[®] Technology



Why have you been given Siltac or Silvex ostomy seals?

These new silicone based ostomy seals may have been recommended to you for a number of reasons:

- You are experiencing leakage
- You have mild, moderate or severe skin problems around your stoma (even slightly red skin can still be a skin problem)
- You are having problems getting a perfect fit of your bag or moulding your current seal
- You are having problems removing the residue left on your skin when changing

What are Siltac and Silvex silicone ostomy seals?

Siltac and Silvex are ostomy seals made from a revolutionary new patented, Responsive Silicone® material that protects your skin completely, allows it to breathe, and leaves no residue or mess.

Why are they different from other ostomy products?

Current ostomy products are made from hydrocolloid. Hydrocolloid seals or rings absorb liquid (and body waste!) to keep it away from the skin.

Level of absorption therefore becomes one of the most important features of a hydrocolloid. This absorption of moisture (and body waste) is important for your current product to stick to wet skin.

We believe that having an adhesive containing faeces and urine on your skin for lengthy periods of time is not a hygienic or effective solution to protect against leakage or keep your skin healthy.

Trio Siltac® and Trio Silvex® are soft silicone ostomy seals that work completely differently. They do not absorb moisture or body waste; for the first time you can have a clean and hygienic seal.

Siltac and Silvex therefore **do not erode, do not leave any residue** and come off **cleanly** and **pain free**.



Comparing Siltac and Silvex to your current product

Siltac and Silvex

Current rings and bags

Always a perfect fit



No Residue



Painless removal



Moves with your stoma



More hygienic



Does not erode/ longer lasting



Instant adhesion - no waiting or warming



Allows your skin to breathe



How does it work?

01

Red/ sore peristomal skin due to contact of output with the skin



02

Siltac/ Silvex hugs your stoma to create a custom fit, always protecting the skin from harmful output



03

By offering complete protection and no scrubbing off residue the skin heals itself and remains healthy




How to use - Important!


Siltac and **Silvex** hug your stoma, so ensuring you have a firm and snug fit.

Do not use Siltac and Silvex in the same way as your current products. Please follow the below instructions.


01 Gently peel the edge of the seal using thumb and forefinger then carefully but firmly pull from the blister pack.

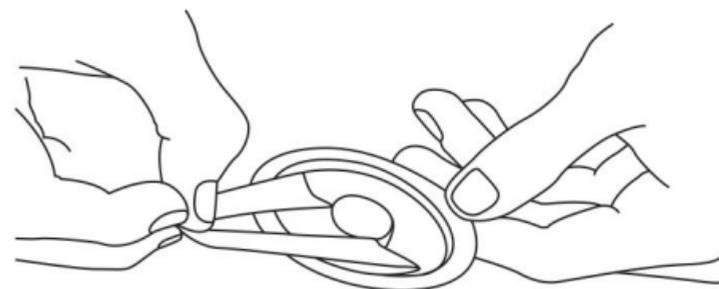

02 Stretch the seal over the stoma and allow it to retract back, forming a perfect seal around the stoma

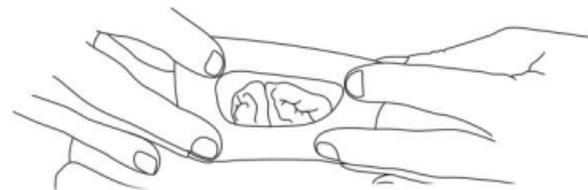

03 Once the seal has retracted around the stoma, gently stick it to the skin

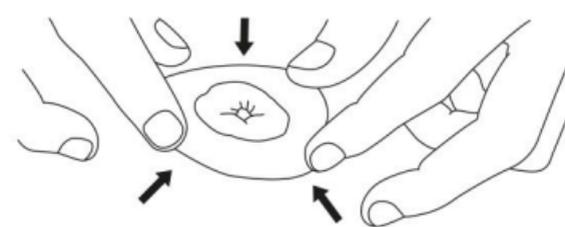

PLEASE NOTE: If you experience any problems fitting the seal try the next size down to achieve a more secure fit or contact the details on the reverse side of this leaflet.

Important: If there are any gaps between the seal and stoma, gently lift that part of the seal and allow it to hug the stoma, removing any gaps.

Available now - Ordering Information

Trio **Siltac**® - Silicone ostomy seal

30 silicone cohesive seals (20mm-28mm) REF TR 1020

30 silicone cohesive seals (28mm-35mm) REF TR 1028

30 silicone cohesive seals (35mm-44mm) REF TR 1035

Trio **Silvex**® - Silicone convex seal

10 silicone convex seals (20mm-30mm) REF 1020CX

10 silicone convex seals (30mm-40mm) REF 1030CX

Trio **Sillex**® - Silicone flange extender

20 silicone flange extenders REF TR1060

Trio **Silken**® - Silicone Stoma Gel

60gm tube - silicone stoma gel REF TR1070

Trio **Elite**® - Sting free adhesive remover

30 wipes REF TR202

50ml aerosol REF TR201

Trio **Elisse**® - Sting free skin barrier

30 wipes REF TR203

50ml aerosol REF TR204

Trio **Pearls**® - Super absorbent - gelling & odour control sachets

Pearls REF TR205

Also available at www.TrioOstomyCare.us





Follow us on  @TrioOstomyCare  Trio Ostomy Care

Professional Products & Services LLC
USA Distributor
PO Box 818 Loughman, FL 33858
1-877-400-1700 · 863-421-9410 (Fax)

E: mail@trioostomycare.us

TrioOstomyCare.us